

DL8000 bp DNA Ladder

Code No. :	NG109S	250 µl×5 (250 T)
	NG109M	250 µl×10 (500 T)

Introduction :

This product is a ready-to-use DNA molecular weight marker pre-mixed with 1x loading buffer, consisting of 6 linear double-stranded DNA. Suitable for the estimation of size and rough quantification of double-stranded DNA sizes ranging from 100 bp to 8 kb.

The 6 bands of this product are 100, 250, 500, 750, 1000, 2000, 3000, 5000 and 8000 bp, with concentrations of 100 ng of 750 bp, the other is 50 ng, when loading 5 µl.

Storage buffer:

10 mM TrisCl (pH 8.4) , 10 mM EDTA , 0.02 % bromophenol blue , 5 % glycerol.

Storage :


Store at -20 °C for up to a years and avoid repeated freezing and thawing. RT for 2-3 months.

Usage:

1. Recommended for agarose gel electrophoresis at 1.0~2.0 %; not recommended for polyacrylamide gel electrophoresis;
2. Electrophoresis buffer can be 1×TAE or 0.5~1×TBE; voltage 6~8 V/cm gel length, electrophoresis time 30~60 minutes; voltage 20~30 V/cm gel length, electrophoresis time 10~15 minutes;
3. Using a sterile pipette tip, aspirate 5~10 µl of this product and add it to the loading well according to the wells;
4. Start electrophoresis after adding the DNA sample to be tested;
5. After electrophoresis, stain with ethidium bromide (EB) or another DNA dyes and observe the bands.

Precautions :

1. Testing indicates that the band pattern remains unchanged when stored at room temperature for one month; however, low-temperature storage is recommended to prevent band degradation used by nuclease contamination due to improper handling.
2. Do not heat before use.
3. Replace the electrophoresis buffer promptly when its buffering capacity decreases to avoid affecting the resolution.

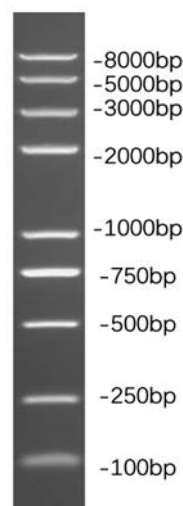

Load 5 µl onto a 1% TAE agarose gel, EB staining